

Fisherbrand™ Mini-Centrifuge

Mini but Mighty - ideal for small runs

The Mini Centrifuge is ideal for microfilter cell separations and HPLC samples, with adapters to accommodate all major microtube and PCR tube styles (1.5/2.0, 0.5, 0.2 mL) for customised use.

Easy to Use

Tool-free SnapSpin rotors and quick-spin options make the centrifuges fast and easy to use.

Ergonomic Design

- The lid fits in the palm of the hand for one-handed locking, unlocking, and quick tube removal. The lid height accommodates taller lidded and skirted tubes
- The ABS plastic housing and polycarbonate plastic lid make cleaning easy

Safety Measures

- Fitted with a braking function, the unit quickly stops when powered off or when lid is opened
- A weighted rubberised base eliminates potential movement and dampens vibration and noise



Accessories and Storage Case Included


Specifications:

Maximum Speed: Approx. 6,000 rpm (fully loaded)
Maximum RCF: 2,000 x g
Power: 100–240 VAC 50/60 Hz 0.6A, Approx. 7.3 W
Five-year warranty



Package includes:

- Mini-Centrifuge
- 6-Place 1.5/2.0 mL Tube Rotor
- 16-Place 0.2 mL Tube Rotor (singles, or 2 x 8 strip)
- 6 Tube Adapters 0.2 mL
- 6 Tube Adapters 0.5 mL
- 4 O-Rings
- Storage Case for Rotors and Adapters
- Tube Rack
- Low Voltage, Double-Insulated Power Adapter with 4 Interchangeable Plugs

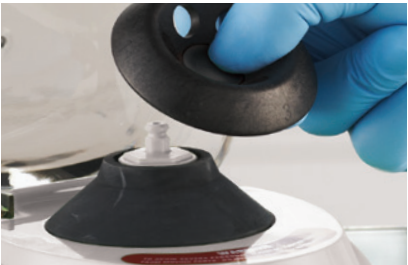

SnapSpin Rotors change out with the push of a button. No tools are needed

Cat. No.	Colour	L x W x H cm	Qty
16617645	Gray	15.3 x 12.8 x 10.9	1

Replacement Parts	Qty
12883332 SnapSpin Rotor 6 x 1.5/2.0 mL Tube	1
12893332 Snapspin Rotor 2 x 8 Strip 0.2 mL	1